

QC now

THE BEST WAY TO SECURE
INTERNAL QUALITY CONTROL



THE INTERNAL QUALITY CONTROL SOLUTION FOR YOUR DIATRON PLATFORMS

SOLVE QC PROBLEMS EFFICIENTLY

AND CONCENTRATE ON TESTS THAT REALLY
REQUIRE YOUR ATTENTION



SAVE TIME
by avoiding
false positives

QC now is designed to meet your needs. The QC add-on provides long-term statistical analysis of internal control results coming from the control software of STRATEC/Diatron devices such as the P780 clinical chemistry analyzer or the KleeYa chemiluminescence immunoassay analyzer.



SAVE COST
reduce time delays
by providing real
time data analysis



HIGH SERVICE QUALITY
full daily statistics at a glance
and secure QC records using
an external database

FEATURES AND BENEFITS

AUTOMATIC OUTLIER DETECTION
see the tests that require your attention

**VIOLATION DIALOG
FACILITATING COMPLIANCE**
to review rule violations and document corrective actions

AUTOMATE QC DATA TRANSFER
utilize real-time access to QC data to improve trouble-
shooting and eliminate delays in QC data submission deadlines

AUDIT TRAIL WITH FILTER FUNCTION
for total traceability of user interaction


CUSTOMIZED ANALYSIS DECISION CRITERIA
select the decision criteria to make the analysis of the data
point

OPTIMIZED DATA INTERPRETATION
easy-to-read charts and tabular views that support you in
data interpretation and decision making



QC REPORT AT ANY TIME
monitor and improve the effectiveness of your QC and
review reports online or print them locally for efficient
QC management

**REMOTE CONNECTIVITY
AND "SILENT INSTALL"**
provides simple installation without users having to
interact

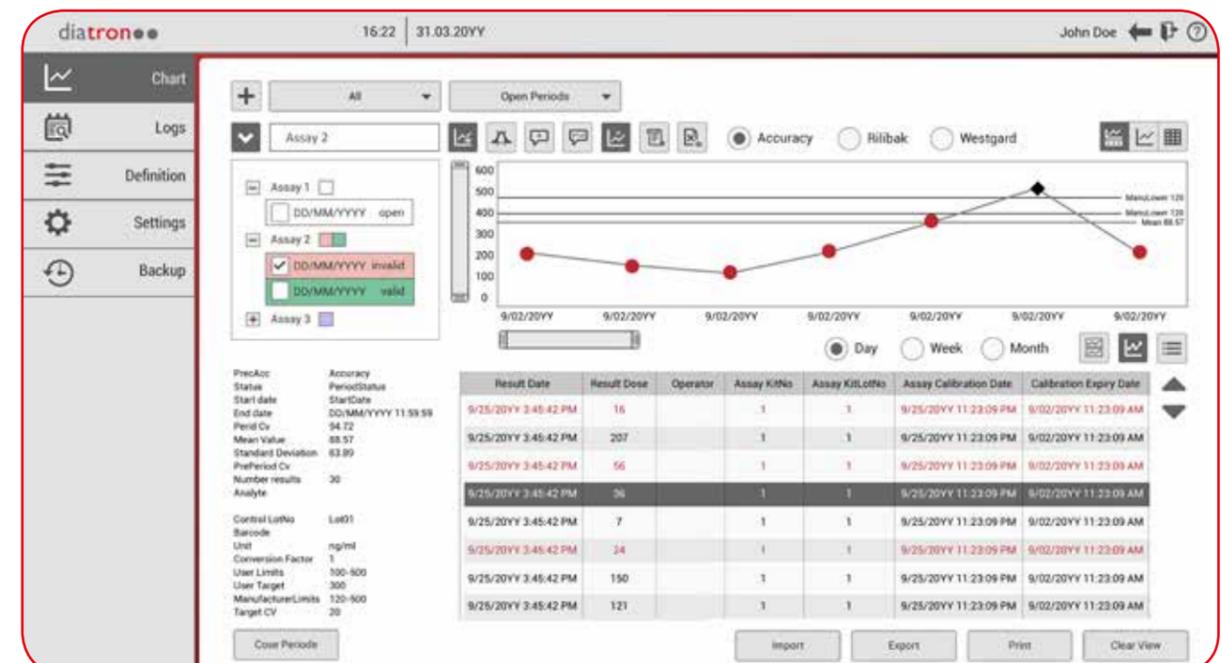

SPECIFICATIONS

CHARACTERISTICS

Application	Control result evaluation based on target value range or “intra laboratory” limits
Access rights	Two user groups: Administrator and Operator
Compatibility	KleeYa, P780 Other STRATEC / Diatron devices on request
Control types	Westgard RiLiBÄK Accuracy Precision
QC Trending	Levey-Jennings daily and day-to-day plots Youden plots
Report formats	PDF CSV Hard copy
Report types	
<i>QC Period Report</i>	Displays all information related to the data generated with the QC software
<i>Audit Trail Report</i>	Displays all events leading to the creation, modification and deletion of an electronic record
<i>Event Report</i>	Displays all messages from the system to the user, in particular error messages
Data transfer	Automatic collection of data directly from analyzer main software or manual import via QC software UI
Archive	Remove closed periods from the database and move them to a configurable path
Backup	Create local compressed and stored copies of files and folders
Regulatory	Compliant with RiLiBÄK



We are a member of the STRATEC Group | www.diatron.com | sales@diatron.com

Diatron MI Zrt makes best efforts to provide accurate information and specifications herein, whereby it may be possible for omissions or errors to have occurred. Diatron MI Zrt shall have no liability for the accuracy and completeness of the information and cannot be held liable for any representations or warranties, expressed or implied, and any third party claims or losses of any damages. Diatron MI Zrt reserves the right to make any changes to the information, including specifications, without notice. The applications and products herein may not be available in all markets. ©2019, Diatron MI Zrt, Hungary, all rights reserved.

Version Nr QC200520VI