


diatron●●

P500

INTELLIGENT WALK AWAY
CHEMISTRY

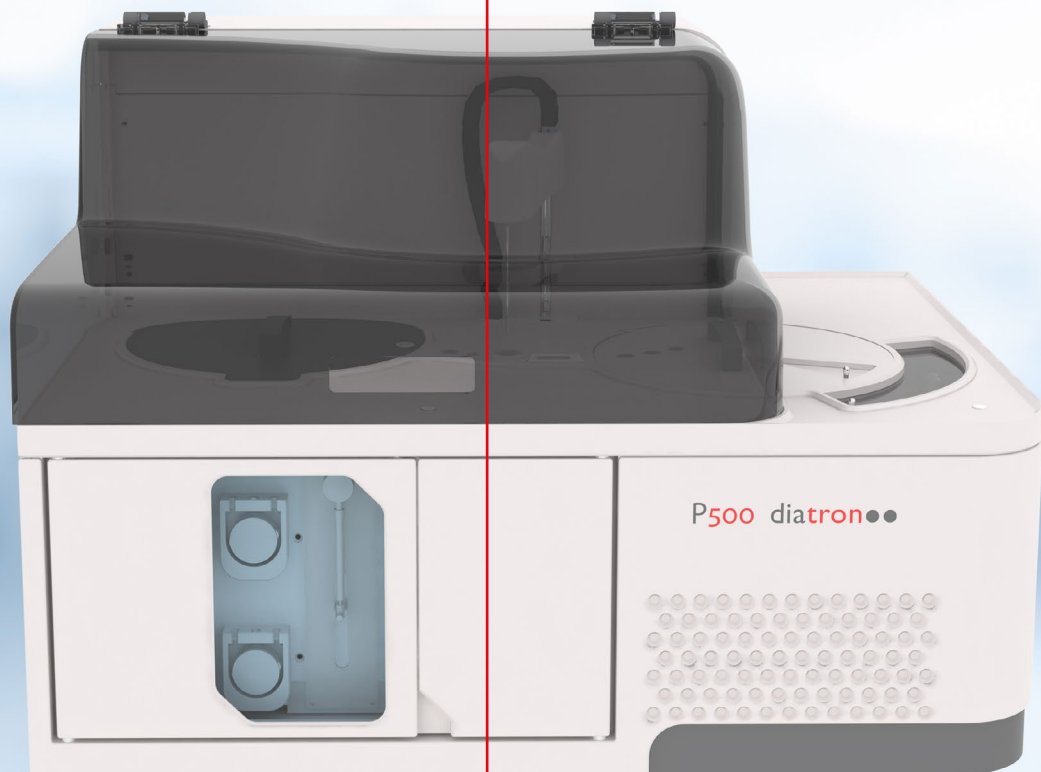

INTELLIGENT WALK AWAY CHEMISTRY

P500

The Diatron P500 is an ergonomically designed, fully automated, random access bench top clinical chemistry system. As a result of its innovative design and advanced Windows-based software, the P500 can achieve a total throughput of up

to 420 tests per hour. The P500 offers features usually only higher cost systems have such as uninterrupted workflow, enhanced walk-away operation and remote access diagnostics, which positioning it as a market leading system.

FLEXIBILITY

State of the art benchtop system with removable reagent tray, host-query LIMS, remote access capability & diagnostics and intuitive, windows based software

PRODUCTIVITY

Continuous sample and reagent loading with on board capacity of 95 samples (5 racks x 19 position) and up to 72 reagent vials + ISE

RELIABILITY & TRUST

Diatron manufactures high performance analysers with best available quality components



SPECIFICATIONS

TECHNICAL SPECIFICATIONS

- Electrical requirements: 600 VA
- Weight: 115 kg, 253 lbs
- Dimensions (Width x Depth x Height): 90 cm x 66 cm x 62 cm
35" x 26" x 24"
- Water Consumption (l/h): ~ 1.5l
(up to 2l depending on usage)
- Computer/SW environment: up to Win 10
- Interferential optical filters (340-750 nm): 12
- Photometric range (linearity): -0.1/3.6
- Optical resolution: 0.0001 ABS
- Optional direct 3-ISE module
(installable by service)
- Internal Barcode Reader for samples, sectors and reagents
- Halogen Tungsten lamp life with sleep mode 2000 h
- Sample and reagent precision of syringe at full stroke: 0.05 % CV

GENERAL CHARACTERISTICS

- Small to mid-size labs random access analyser (fully automated)
- Maximum throughput with ISE: 420
- Temp. controlled cuvette tray via air bath (min maintenance): 37 ± 0.1 °C
- Reaction volume: 180-650 µl
- Sample volume: 2-200 µl
- Cuvette Path/no. of cuvettes on board: 6 mm/80 units
- Probes/Pipettes: 1
- Time to first result (seconds): 45

GENERAL CHARACTERISTICS

- Removable and interchangeable reagent tray
- Refrigerated reagents: 8 ± 2 °C
- Refrigerated on-board single reagent vials: 72
- Up to three different reagents per test
- Continuously monitored loading/access of reagents
- Multiple reagent bottles of the same assay/test can be on-board (different lots)
- On-screen inventory management with real time tracking of on-board reagent availability
- Software available in 9 languages

SAMPLE HANDLING

- Primary, secondary tubes and micro-cups for pediatric samples
- Initial rack quantity: 5 x 19 samples with continuous loading up to 99 racks providing a capacity of 1100 samples per day
- Samples on-board capacity: 95
- Continuous loading of sample racks
- Real STAT function
- Automatic sample dilution
- Post-analytical tasks such as re-running and post-dilution
- Editable complete A4 patient report
- Unlimited/searchable patient records

Analyzer and ISE module

SAP No	Description	Catalogue No
10025007	P500 Analyzer without ISE	P500-1
10025008	P500 Analyzer with direct 4-ISE module	P500-2
10021291	Pictus ISE module (module for upgrading Pictus and P models)	5008
10030037	P500 Analyzer without ISE (open)	P500-1-O
10025009	P500 Analyzer with direct ISE module (open)	P500-2-O

diatron●●

We are a member of the STRATEC Group | www.diatron.com | sales@diatron.com

Diatron MI Zrt makes best efforts to provide accurate information and specifications herein, whereby it may be possible for omissions or errors to have occurred. Diatron MI Zrt shall have no liability for the accuracy and completeness of the information and cannot be held liable for any representations or warranties, expressed or implied, and any third party claims or losses of any damages. Diatron MI Zrt reserves the right to make any changes to the information, including specifications, without notice. The applications and products herein may not be available in all markets. ©2025, Diatron MI Zrt, Hungary, all rights reserved.

Version Nr. P500250615USV7

